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
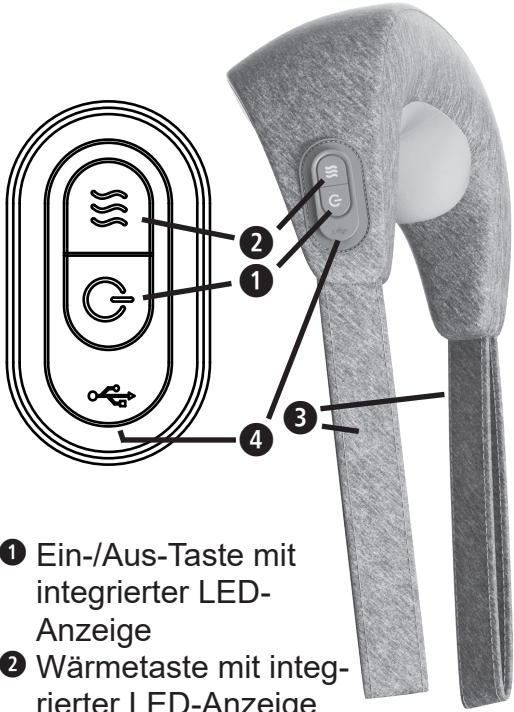
Neck Massager NM 900

Device and controls

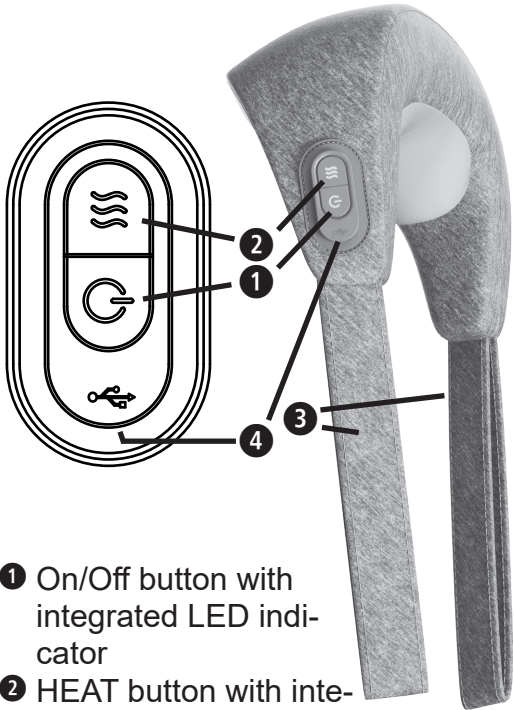- ❶ On/Off button with integrated LED indicator
- ❷ HEAT button with integrated LED indicator
- ❸ Hand strap (each 1 on right and 1 on left side)
- ❹ Charging port (USB connector)

Explanation of symbols


Do not puncture the pad!


Not suitable for children under 3 years!

The device must only be used indoors!

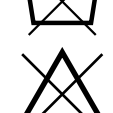Do not wash!

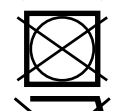Do not bleach!

Do not dry the article in a tumble drier!

Do not iron!

Do not dry clean!

IMPORTANT
Non-observance of these instructions can result in serious injury or damage to the device.

WARNING
These warning notes must be observed to prevent any injury to the user.

CAUTION
These notes must be observed to prevent any damage to the device.

NOTE
These notes give you useful additional information on the installation or operation.

LOT
LOT number

Manufacturer

Recycling symbols / codes:
These are used to provide information about the material and its proper use and recycling.

GB

Important Instructions

Retain for future use

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

- The device must only be used for its intended purpose as described in the instruction manual. Using the device for any other purpose invalidates the warranty.

- The device is not to be used for commercial use or medical purposes.
- If you have health concerns, consult your doctor before using the device.

Only use the device indoors.

Do not use the device in damp rooms (e.g. when bathing or showering).

- Keep the device away from direct sunlight, do not expose it to high temperatures.
- Follow all legal constraints and regulations for disposing of the device and the internal battery. (See also Section “Disposal”.)

Safety Instructions – About the power supply

- This device must only be supplied at safety extra low voltage (SELV) corresponding to the marking on the appliance.
- Do not carry, hold or drag the device by the cable!
- Batteries must never be squeezed, disassembled or affected by hot temperatures (above 60°C). There is the danger of fire, explosions and the risk of a leakage of dangerous and poisonous chemicals.
- This device contains non-replaceable batteries.
- The batteries contained in this device may be replaced only by service technicians that have been authorised by **medisana**.
- The battery life gets reduced if it is not used for a long time. It is therefore recommended to charge the battery at least once or twice per month.
- Use a USB mains adapter with an output of 5 V and (at least) 2 A for charging the battery. A larger output voltage will damage the device.
- Only insert the power supply unit into a socket when the device is switched off. The device must not be used while it is being charged.
- Do not touch the mains adapter if you are standing in water or if your hands are wet or damp. Otherwise, there might be the risk of damages or an electric shock.

Safety Instructions – Specific groups of people

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children must be supervised to ensure that they do not play with the device.
- Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.

- Do not use the **NM 900** neck massager, or only after consulting your doctor, if:
 - »you are pregnant,
 - »you have a cardiac pacemaker, artificial joints or electronic implants,
 - »you suffer from one or more of the following illnesses and complaints: Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- Do not use the appliance near your eyes or near any other sensitive areas.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- The surface of the appliance is hot. People who are insensitive to heat must be careful when using the appliance.
- Before using the massager, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.
- Children under the age of three must not use this device due to their inability to react to overheating.

Safety Instructions – Using the unit

- Always carefully inspect the device for damage before use. Do not attempt to use a defective device.
- Do not use the device if it shows signs of damage, if it is not in perfect working order, or if it has been dropped or become damp. In order to prevent hazards, always send the unit in to the service centre.
- Only use the device in the proper place, that is, around the neck and on the shoulders.
- Do not use the device while you are asleep or when you are lying in bed.
- Do not use the device before going to sleep. The massage has a stimulating effect.
- Never use the device while driving a vehicle or operating a machine.
- The maximum operating time for one session is 15 minutes. Prolonged use reduces the service life of the device.
- Very long periods of use can lead to overheating. Always allow the device to cool down before using it again.
- Never put the unit down next to an electric stove or other sources of heat.
- Do not allow the massage cushion to come into contact with pointed or sharp objects.
- Do not cover up the appliance when it is switched on. Never use it under blankets or cushions. There is a risk of fire, electric shock or injury.

Safety Instructions – Maintenance and cleaning

- The device is maintenance free. Do not try to open it. There are no parts inside that are to be cleaned or serviced by the user.
- You may only clean the device yourself. If the device malfunctions please do not repair it yourself as this leads to the automatic expiry of all warranty rights. Consult your specialist dealer and have repairs carried out by the service centre to avoid dangers.
- Do not immerse the device in water or other fluids. If, nonetheless, liquid should penetrate the device, unplug the mains cable immediately.
- Cleaning and user maintenance must not be carried out by children without supervision.

Scope of delivery

Please check first whether the unit is complete and does not show any damage.

If in doubt, do not operate the unit and contact your dealer or service centre.

The scope of delivery includes:

- 1 **medisana** Neck Massager **NM 900**
- 1 Charging cable, USB-C to USB-C (cable length: approx. 1 m)
- 1 Quick Start Guide (short instructions)

If you notice any transport damage when unpacking, please contact your dealer immediately.

Packaging is recyclable or can be returned to the raw material cycle. Please dispose of packaging material that is no longer required properly.


WARNING
Please ensure that the packaging material is kept away from the reach of children! Risk of suffocation!

Charging the battery

- Use a 5 V mains adapter with an output current of (at least) 2 A. Connect the supplied USB charging cable to the adapter and to the USB port ❹ of the device.
- The device must not be used during the charging process.
- The On/Off button ❶ indicates the charging process by flashing in red colour. When the battery is fully charged, the button ❶ is lit continuously in green colour.
- If the On/Off button ❶ is flashing in red colour the remaining charge is low – the battery should be recharged as soon as possible.

Using the device

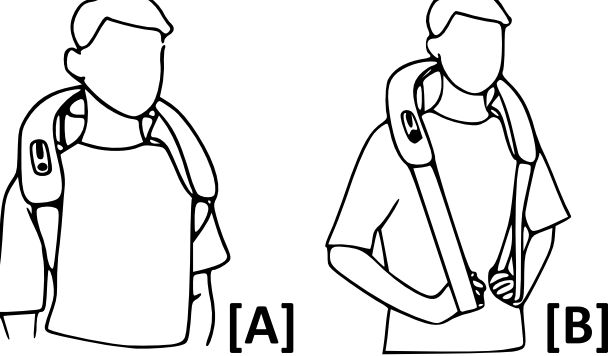Place the massager around the neck **[A]** so that it has sufficient hold on the shoulders.

- Press the On/Off button ❶ for (approx.) 2 seconds to switch the device on. It starts immediately with the massage. The operation is indicated by the On/Off button ❶ being lit in green colour.
- A short press on the Heat button ❷ switches the heat function on or off. The button ❷ is lit in red colour when the heat function is active.
- A short press on the On/Off button ❶ changes between two speed levels of the massage heads. Besides that, the device changes the massage direction automatically every 2 minutes.
- Put your hands through the hand straps ❸ on the right and left sides **[B]**. By pulling on the straps you may stabilize the device in its position and adjust the massage intensity (pulling stronger increases the massage intensity).
- Press the On/Off button ❶ again for (approx.) 2 seconds to switch the device off. If you do not switch it off manually it switches off automatically after (approx.) 15 minutes.

With a fully charged battery, the device can be used for approx. 1 hour, depending on the speed level and whether heat is used.

Cleaning and maintenance

- Before cleaning the unit, make sure it is switched off and disconnected from the mains. Remove the charging cable.
- Before cleaning the appliance, make sure that the appliance is switched off and the power supply unit is disconnected. Allow the unit to cool down.
- Clean the massager with a slightly damp sponge only. Never use brushes, aggressive cleaning agents, petrol, thinners or alcohol. Rub the unit dry with a clean dry cloth.
- May not be bleached, ironed or dry cleaned!
- Never immerse the appliance in water to clean it and ensure water cannot get into the appliance.
- Wait until the device is completely dry before using it again.
- You should ideally store the device in its original packaging and keep it in a clean, dry place.
- Store the device in a safe, clean, cool and dry, well-ventilated place. Secure it against falling and shocks.

Technical Data

Name and model:	medisana Neck Massager NM 900
Power supply:	Li-Ion battery 7.2V, 2000mAh input: USB-C, 5V=, 2,0A approx. 2 hours.
Charging time:	
Max. operating time:	15 minutes, with automatic switch-off
Dimensions:	approx. 43.5 × 23 × 12 cm
Weight:	approx. 0.6 kg
Operating conditions:	Only use in dry rooms at –5°C to +40°C
Storage conditions:	In a cool, dry place

Disposal

 This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether  or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

Warranty and repair terms

Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:


- The warranty period for **medisana** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
- Defects in material or workmanship will be removed free of charge within the warranty period.
- Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
- The following is excluded under the warranty:
 - All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - Accessories which are subject to normal wear and tear.
- Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

medisana GmbH, Carl-Schurz-Str. 2, 41460 NEUSS, GERMANY

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

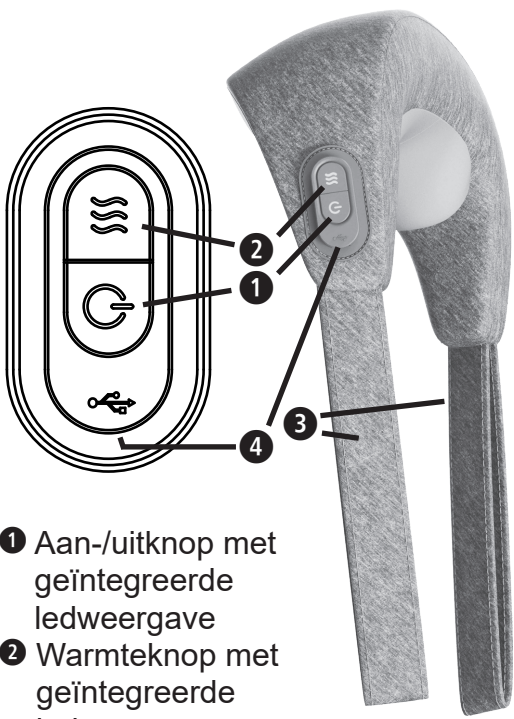The current version of this instruction manual can be found under <https://docs.medisana.com/88965>.
Information about service can be found here:
www.medisana.com/servicepartners

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