
medisana®

DE

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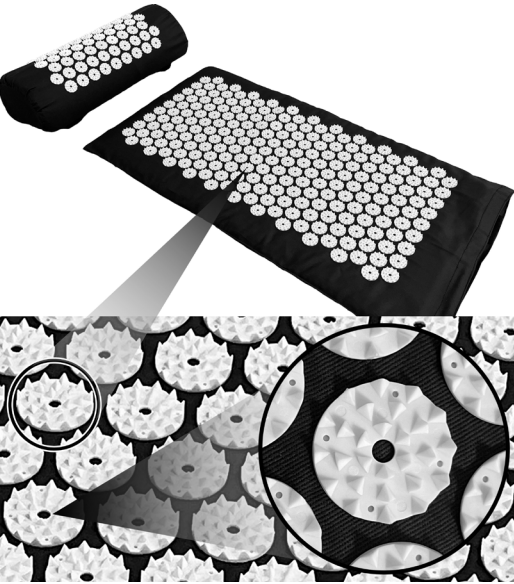
Acupressure mat with heat MM 100

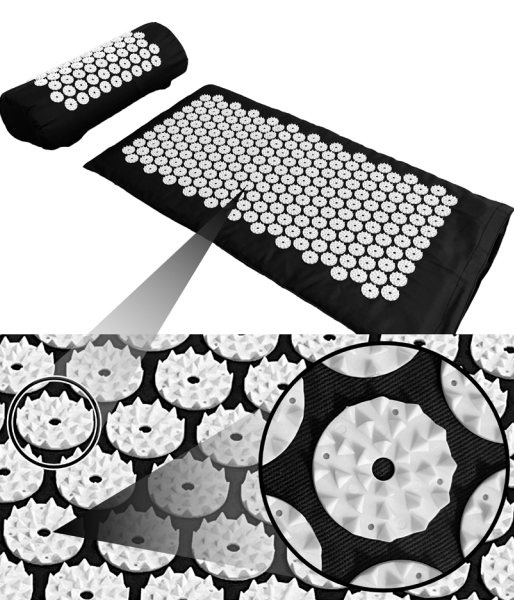Explanation of symbols

Do not use the device when it is folded!

Do not puncture the device!


Not suitable for children under 3 years!


The outer cover can be cleaned in a gentle wash cycle at max. 30°C!

Do not wash the acupressure mat, nor the (inner part of) the pillow!

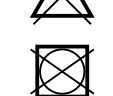Do not bleach!

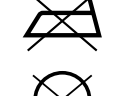Do not dry in a tumble drier!

Do not iron!

Do not dry clean!

The device must only be used indoors!

IMPORTANT
Non-observance of these instructions can result in serious injury or damage to the device.

WARNING
These warning notes must be observed to prevent any injury to the user.

CAUTION
These notes must be observed to prevent any damage to the device.

NOTE
These notes give you useful additional information on the installation or operation.


Protection category II

LOT number

Manufacturer

Recycling symbols / codes:
These are used to provide information about the material and its proper use and recycling.

GB

Important Instructions

Retain for future use

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

- The device must only be used for its intended purpose as described in the instruction manual. Using the appliance for any other purpose invalidates the warranty.
- The device is not to be used for commercial use or medical purposes. It must not be used in hospitals.
- When the device is likely to be used for an extended period, controls should be adjusted to a setting recommended for continuous use.

Safety Instructions – About the power supply

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- The control switch should not be covered or placed on or below the mat while the unit is in operation.
- Keep the mains lead away from hot surfaces.
- Never carry, pull or turn the device by the mains lead and do not let the lead get tangled.
- Never touch a device that has fallen into water. Unplug the unit from the mains outlet immediately.
- The controls and leads must not be exposed to any kind of moisture.
- The device may only be operated with the control unit (SS43A) supplied.

Safety Instructions – Specific groups of people

- This appliance must not be used by persons insensitive to heat and other very vulnerable persons who are unable to react to overheating.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Children below the age of 3 years are not allowed to use this device, because they are not able to react on overheating.
- The device must not be used by young children over the age of 3 years, unless the control unit has been adjusted accordingly by its parents or legal guardians and unless the child has been well instructed how to use the control unit properly.
- The appliance can be used by children older than 3 years and younger than 8 years under supervision and with the control always set to minimum temperature value.
- Should you have any concerns about health issues, consult your doctor before using the device.
- Do not use the device on parts of the body that are swollen, injured or inflamed.
- Should you experience persistent pains in the muscles and joints, please inform your doctor. Persistent pains can be a symptom of serious illness.
- If you wear a pacemaker, please be aware that the electromagnetic fields emanating from this electric product may, under certain circumstances, interfere with the functions of your pacemaker. Therefore, ask your doctor and the manufacturer of the pacemaker for advice before using this product.

- If you experience any pain or discomfort while using the unit, stop using it immediately.
- The cable and control unit of the appliance can give rise to risks of entanglement, strangulation, tripping or treading if not correctly arranged. The user shall make sure that excess ties and electric cords shall be arranged in a safe way.

Safety Instructions – Operating the device

- Check the article carefully for signs of wear and/or damage before each use. Do not operate the device if you notice wear, damage or signs of improper use, or if it does not function correctly. Instead, contact an authorised service centre in such a case.
- Do not use if you notice wear, damage or signs of improper use on the mat, the switch or the cables or if the device is not working.
- The control unit must not be placed on or under the mat during operation and it must not be covered.
- Do not cover the mat with any other cushions.
- Do not fold or crease the acupressure mat during use.
- Neither pierce the mat with, nor attach safety pins or other spiky or sharp objects to the mat.
- Never use the heating function unsupervised.
- Do not fall asleep while the heating is switched on. Prolonged use of the heating function at a high setting could lead to skin burns.
- Ensure that you follow the safety notes, as prolonged use at a high setting can result in a heat stroke or burns to the skin.
- Do not sit on the acupressure mat. Instead, rest the mat on or against the part of the body to be treated. The intended use of the is to treat only back and neck.
- Do not use the device when wet and only use it in dry environments (not in the bathroom or similar environments).

Safety Instructions – Maintenance and cleaning

- You are only permitted to clean the unit. To avoid hazards, never repair the device yourself. Contact the customer service.
- Cleaning and user maintenance shall not be made by children without supervision.
- If the power supply cable is damaged, it may only be replaced by **medisana**, an authorised dealer or by qualified personnel in order to avoid hazards.
- A damaged power supply cable must be replaced with a power cable of the same type.
- If a fault occurs, do not attempt to repair the unit yourself. Repairs must only be carried out by an authorized specialist dealer or other suitably qualified personnel.
- When storing the appliance, allow it to cool down before folding.
- Do not crease the appliance by placing items on top of it during storage.

Scope of delivery

Please check first whether the unit is complete and does not show any damage.

If in doubt, do not operate the unit and contact your dealer or service centre.

The scope of delivery includes:

- 1 **medisana** Acupressure mat with heat **MM 100** with removable control unit **SS43A** and head pillow
- 1 storage pouch
- 1 Quick Start Guide (short instructions)

If you notice any transport damage when unpacking, please contact your dealer immediately.

 Packaging is recyclable or can be returned to the raw material cycle. Please dispose of packaging material that is no longer required properly.


WARNING

Please ensure that the packaging material is kept away from the reach of children! Risk of suffocation!

Control Unit SS43A

- ① LED display of the heating level (available: levels 1...9)
- ② Button to increase heating level
- ③ On/Off button
- ④ Button to reduce heating level

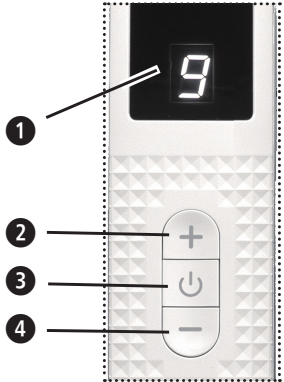Using the device

Make sure, the connection cable of the control unit is properly connected to the plug of the acupressure mat. Plug the unit into the mains outlet and press the On/Off button ③. The device starts at lowest level (level 1), which is displayed at ①. After a few minutes you will be able to feel the mat getting noticeably warmer. If you require a higher temperature setting, press the “+” button ② several times until you reach the required level. (The maximum level is level 9.) If the acupressure mat gets too hot for you, press the “–” button ④ several times until you reach the required level. After a short time, the mat will cool down. The acupressure mat may cool down after a short period of time. This is normal. The fast re-heating according to the setting will follow immediately.

For continuous use, you have to always select the lowest level 1.

To switch the heating off, press the On/Off button ③. The LED display ① will go out, indicating that the unit is switched off. Unplug the unit from the mains outlet.

The heating will switch off automatically after approx. 90 minutes of continuous use. Unplug the unit from the mains outlet once you have finished using it.

Troubleshooting

If there is no function at all after switching on the device (no heating level is displayed at ①), please check if the mains plug is correctly plugged in the mains outlet. If this is the case, the device is defect. Then please do not use the device and contact the service center.

Cleaning and maintenance

- Before cleaning the mat, unplug it from the mains outlet and let it cool down for at least ten minutes.
- The acupressure mat and the head pillow are equipped with a removable outer cover.
- You may wash the outer covers (i. e. without pillow, without mat, without control unit) gently by washing machine. But note that the mat may bleed some colour when washed, and the plastic tiles may rub against other laundry. Therefore, we recommend to wash the covers separately.
- However, the best way is to put each cover in a bath with lukewarm water and some mild detergent and then gently squeeze it out. Rinse the cover several times to remove all the detergent.
- Never use strong detergents or hard brushes.
- To dry, lay the covers flat on a moisture-absorbent towel or mat. Re-fasten the covers to the mat and pillow only once they are completely dry.
- You should untangle the cable from time to time.
- Store the unit in a clean and dry place, laid out flat with nothing covering or resting on top of it.

Non-washable parts:

- Before cleaning the unit, make sure it is switched off and disconnected from the mains.Clean the unit once it has cooled down using a soft, damp cloth.
- Under no circumstances should you use aggressive detergents, strong brushes, solvents or alcohol, as they may damage the surface of the unit.
- When cleaning, never immerse the unit in water and make sure that no water gets into the unit.
- Wait until the device is completely dry before using it again.
- Undo any kinks in the cable.

Technical Data

Name and model:	medisana Acupressure mat with heat MM 100 (with control unit SS43A)
Power supply:	via Control Unit SS43A – input: 220–240V~ 50/60Hz
Heat output:	approx. 100 W (at max. heating level)
Heating levels:	9 levels. The active level (1 ... 9) is displayed via LED display ①
Max. operating time:	approx. 90 minutes, with automatic switch-off
Dimensions (L×W×H):	approx. 74 × 40 × 2 cm
Weight:	approx. 0.7 kg
Operating conditions:	Use only in dry indoor areas
Storage conditions:	Laid out flat and dry
Item number:	88977
EAN number:	40 15588 88977 6

Disposal

 This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

Warranty and repair terms

Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is. The following warranty terms apply:

- 1.The warranty period for **medisana** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
- 2.Defects in material or workmanship will be removed free of charge within the warranty period.
- 3.Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
- 4.The following is excluded under the warranty:
 - a.All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b.All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d.Accessories which are subject to normal wear and tear.
- 5.Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

medisana GmbH,
Carl-Schurz-Str. 2, 41460 NEUSS, GERMANY

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

The current version of this instruction manual can be found under <https://docs.medisana.com/88977>.
Information about service can be found here:
www.medisana.com/servicepartners

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