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Device and controls

Explanation of symbols

 **IMPORTANT**
Non-observance of these instructions can result in serious injury or damage to the device.

 **WARNING**
These warning notes must be observed to prevent any injury to the user.

 **CAUTION**
These notes must be observed to prevent any damage to the device.

 **NOTE**
These notes give you useful additional information on the installation or operation.

 Only use the device indoors.

 Do not use the device in damp rooms (e.g. when bathing or showering).

LOT LOT number

 Manufacturer

 Recycling symbols / codes:
These are used to provide information about the material and its proper use and recycling.

GB Safety instructions

 **Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.**

- Only use the device for the according the instructions for use. The guarantee is invalidated if the device is misused.
- The device is not intended for commercial use or medical applications.
- Check the device for damage before use. If in doubt, do not use it and contact customer service.

Power supply

- Before connecting the device to your power supply (for charging), make sure that the mains voltage indicated on the type plate matches that of your mains supply.
- Only charge the device when it is switched off.
- This appliance contains batteries that are non-replaceable.
- The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- Do not reach for a device that has fallen into water.
- Do not carry, pull or twist the device by the charging cable.
- If the cable or device is damaged, the device must not be used. For safety reasons, parts may only be replaced by an authorised service centre. To avoid hazards, send the device to the service centre for repair.
- Make sure that cables do not become a trip hazard. They must not be kinked, pinched or twisted.

Special people

- Do not use the device on children, disabled persons, sleeping persons or persons who are insensitive to heat and cannot react to overheating.
- This device can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they are supervised and have been instructed in the safe use of the device and understand the hazards arising from its use.
- Children should be supervised to ensure that they do not play with the device.
- Do not use this device to support or replace medical applications. Chronic conditions and symptoms may worsen.
- You should refrain from using the massage device or consult your doctor first if:
 - you are pregnant,
 - you have a pacemaker, artificial joints or electronic implants.
 - you suffer from one or more of the following conditions: circulatory disorders, varicose veins, open wounds, bruises, skin tears, phlebitis.
- Use the device exclusively as described in these instructions, with the hips positioned for massage.
- If you experience pain or find the massage uncomfortable, stop using the device and consult your doctor.
- If you have unexplained pain, are under medical treatment and/or use medical devices, consult your doctor before using the massage device.
- Do not use the device while sleeping or lying in bed.
- Never use the device while driving a vehicle or operating machinery.
- Never place or use the device directly next to an electric heater or other heat sources.
- Do not leave the device unattended when it is connected to the mains (for charging).
- Avoid contact between the device and pointed or sharp objects.
- Do not cover the device when it is switched on. Never use it under blankets or pillows. There is a risk of fire, electric shock and injury.

Use

- The maximum operating time for one application is 15 minutes. Prolonged use will shorten the service life of the device and may cause muscle tension instead of relaxation.
- Prolonged use at high settings may cause skin burns. Check the skin areas below the application area regularly to avoid burns or blisters.
- Prolonged use may cause overheating. Always allow the device to cool down before using it again.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

Maintenance and cleaning

- You are only permitted to clean the device yourself. In the event of malfunctions or damage, do not repair the device yourself, as this will void any warranty claims. Do not immerse the device in water or other liquids. If liquid does get into the device, switch it off and do not use it until it is completely dry.
- This appliance contains batteries that are only replaceable by skilled persons.
- Cleaning and user maintenance must not be carried out by children without supervision.

Intended use

The device is intended for self-massage of the abdominal region – for relaxation and well-being. It is not intended for medical or therapeutic use. It is not suitable for use in clinics, massage practices, wellness studios, etc. The massage device is not a medical device and is not intended for pain treatment. If in doubt or unsure, consult your doctor before use.

Device and controls

- ❶ Mode key
- ❷ On/Off key
- ❸ Warmth key
- ❹ Elastic band for fastening around the hips
- ❺ USB port for charging via charging cable
- ❻ LED indicator
- ❼ Massage nubs with fabric cover

Scope of supply

First, check that the device is complete and undamaged. If in doubt, do not use the device and contact your dealer or service centre.

The scope of delivery includes:

- 1 **medisana** Belly Massager **BM 100**
- 1 USB charging cable
- 1 Belt
- 1 Quick start guide

If you notice any transport damage when unpacking, please contact your dealer immediately.

 **WARNING**
Keep the packaging film out of the reach of children!
There is a risk of suffocation!

Charge device

Charge the device before use by connecting the supplied USB charging cable to a suitable power source and the device connection ❺. The LED indicator ❻ shows the current status:
Flashing red: Low battery charge, please charge immediately
Slowly flashing red: Battery is charging
Solid red: Battery is fully charged.

Use


Place the belly massager around your hips so that the massage nubs ❷ are positioned on your abdomen. Use the elastic strap ❹ to adjust the device to a comfortable width so that it fits snugly but does not constrict.

- Switch on the device by pressing and holding the on/off button ❷ for approx. 1.5 seconds. The device is activated in standard mode (massage against flatulence and low heat).
- By repeatedly pressing the heat button ❸, you can switch between low heat, high heat and off.
- By repeatedly pressing the mode button ❶, you can switch between massage modes mode 1 'massage for flatulence' -> mode 2 'massage for constipation' and mode 3 'massage for menstruation'.
- To switch off the device, press and hold the on/off button ❷ for approx. 1.5 seconds.


Removing the fabric cover

The fabric cover of the massage nubs ❷ can be removed for cleaning by hand. To do this, proceed as follows:

- Switch off the device.
- Pull the retaining ring upwards:

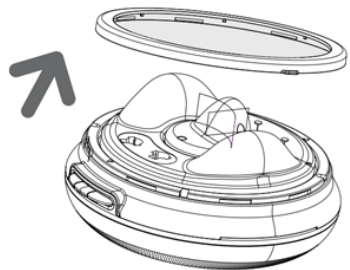Now carefully remove the fabric cover from the retaining ring:

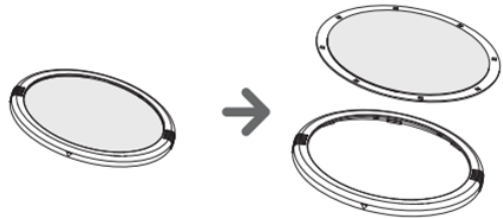Cleaning and care

- Before cleaning the device, make sure that it is switched off and not connected to the power supply (for charging).
- Clean the cooled device only with a damp, soft cloth.
- Never use aggressive cleaning agents, strong brushes, solvents or alcohol. This could damage the surface.
- Never immerse the device in water for cleaning and make sure that no water gets inside the device.
- To clean the fabric cover: Remove it from the device as described, wash it with cold water, allow it to dry and then reinsert it into the device.
- Do not use the device again until it is completely dry.
- It is best to store the device in its original packaging and keep it in a safe, clean, cool and dry place.

Disposal

 This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

Technical specifications

Name and model:	medisana Belly Massager BM 100
Power supply:	7.3V= 2600mAh; Input: 5V= 2.0A
Maximum usage time:	15 minutes
Dimensions L x W x H:	approx. 205 x 205 x 105 mm
Weight:	approx. 1 kg
Storage conditions:	in a cool, dry place
Article number:	88405
EAN number:	40 15588 88405 4

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

The latest version of these instructions for use can be found at <https://docs.medisana.com/88405>.

Service information is available here:
www.medisana.com/servicepartners

Warranty and repair terms

Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **medisana** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. nonobservance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
 5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

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