medisana®

GB Instruction manual Reusable Hot-Cold Pack GP 100 (Head):

IMPORTANT!

This operating manual belongs to this item. The operating manual contains important information regarding use. Failure to follow these instructions may result in serious injury or damage to the product.

Explanation of symbols

CAUTION

These instructions must be followed to prevent possible injury to the product.

This symbol indicates that the product that is normally provided sterile in the same or similar packaging has not been sterilized.

Medical device

Unique product identification

LOT number

Model number: **ACMZ15**

Expiry date

Importer

Manufacturer

Date of manufacture

Authorised representative for the EU

Storage temperature

Recycling symbols/codes: These are used to provide information about the material and its proper use and recycling.

Manual

Service

IMPORTANT INFORMATION! KEEP IN A SAFE PLACE!

Please read this instruction manual carefully, in particular the safety instructions, before using the device. Keep the instruction manual in a safe place for later reference. If you pass the device on to a third party, this instruction manual must remain with the device.

WARNINGS AND PRECAUTIONS

- Only use the product for its intended purpose as described in the instructions for use. Misuse will invalidate the warranty.
- The product is not intended for commercial use.
- Check the product for damage before use – if in doubt, do not use it and contact customer service.
- Always check the temperature before use to avoid burns or frostbite.
- Be careful when removing the pack from the microwave, as it may be very hot. If it is too hot, set it aside until it has cooled down sufficiently to be comfortable to use. Always test the temperature of the pack before each use.
- To avoid frostbite, which can occur if the pack is too cold for your skin, set it aside until it reaches a temperature that your skin can tolerate. If the pack is not cold enough, place it in the freezer for another 20 minutes.
- Do not apply heat/cold to the affected area for more than 15 minutes at a time and take a 20-minute break between applications.
- Do not overheat. Do not leave unattended while heating. Please note that heating for longer than recommended may cause injury or damage to your microwave. Please note that heating times vary in microwave ovens and you may need to adjust the heating time to your preferred temperature.
- Do not heat by any means other than in a microwave.
- Before reheating, check all parts of the pack to ensure that no part is too hot. Knead the pack to ensure that the temperature is evenly distributed.
- Do not wear the product while sleeping. Check the skin regularly during use. If you notice any irritation or burns, remove the product immediately.
- In the case of acute and chronic inflammatory processes, consult a doctor before using the heat compress.
- Do NOT use the product if the contents of the heat/cold pack are leaking and/or the packaging is damaged or torn. If the contents of the heat/cold pack come into contact with your eyes, remove the packaging, rinse the affected area with water and seek medical attention immediately.
- The product should not be used after the expiry date.
- Do not use together with pain-relieving creams, medicinal lotions, creams or ointments.
- People with diabetes, circulatory disorders, nerve damage, paralysis or limited temperature sensitivity should only use the product as directed by a doctor.
- Do not use on sensitive areas and on unhealthy, damaged or broken skin, or on areas where bruising or swelling has occurred within 48 hours, or on body parts that are unable to feel heat or chill or other forms of heat.
- Keep out of reach of children.
- Do not use on individuals who cannot remove the product themselves, such as infants and some elderly people, or who are unable to follow all instructions for use. Children should be supervised by an adult.
- Do NOT apply additional pressure to the product.
- Do NOT overheat the product.
- This product is intended for external use only. DO NOT swallow the product or its contents. Keep out of reach of children and pets.
- Discontinue use and consult a doctor if pain worsens or persists for more than 7 days, or if you experience discomfort, burning, swelling, rash or other changes to your skin.

Intended use

Used to help relieve fevers, bruises, swelling, muscle spasms, aches and pains, headaches, first-degree burns and fatigue. For cold therapy: pain and itching, pre/post exercise, joint pain and sports injuries, swelling and sprains, fever, bruises or swelling, sprains, first-degree burns, after-sun care, headaches and migraines, nose bleeds, relief of fatigue, eyestrain and fatigue, muscle pain and stiffness. For hot therapy: Fatigue Relief, Headache & Migraine, Weary & Puffy Eyes, Muscle Aches & Stiffness, Tired Eyes, Stomachache or Cramps, Aches, Pain & Itching, Pre/Post Workout, Joint Pain & Sport Injuries.

CONTRAINDICATIONS

- People who are allergic to the external material of this product and do not use it on damaged skin;
- Avoid use for those with terminal circulatory dysfunction and those without autonomy;
- Ask a doctor before use if you have diabetes, poor circulation or heart disease, rheumatoid arthritis or are pregnant.

Scope of delivery

First, please check that the product is not damaged in any way. If in doubt, contact your dealer or service centre.

The scope of delivery includes:

- 1 **medisana** Hot Cold Pack **GP 100**
- 1 Instruction manual

If you notice any transport damage when unpacking, please contact your dealer immediately.

CAUTION

Keep the packaging film out of the reach of children!
There is a risk of suffocation!

Cooling / Use as a cold pack:

Place in the freezer for approximately 2 hours before use.

Heating / Use as a hot pack:

- The pack should be at room temperature before heating.
- Place the pack in the centre of a clean microwave (on a turntable if available), without touching the sides. Ensure that all additional microwave functions such as "browning" etc. are switched off.
- Heat the product in the microwave for the time specified below.
- Check that the temperature is suitable before use.
- If you prefer a higher temperature, heat the product for the additional time specified below until the desired temperature is reached. Before reheating, check all parts of the product to ensure that no part is too hot. Knead the product to ensure that the temperature is evenly distributed. Turn the pack 180 degrees before placing it back in the microwave.

IMPORTANT! DO NOT OVERHEAT!

Always reheat from room temperature!

Microwave Rating	Maximum Time
800 W	35 seconds
1000 W	30 seconds

Possible additional time: **10 seconds**

Care instructions

To clean, wipe with a damp cloth. Place on a flat surface and allow to air dry.

Do not wash!

Do not wring out!

Do not tumble dry!

Do not dry clean!

Do not iron!

Dry laid out flat!

Serious incidents must be reported to the manufacturer and the competent authority of the Member State in which the user is established (in Germany, the BfArM).

Technical specifications

Name and model:	medisana Hot Cold Pack GP 100 , model: ACMZ15, for a head circumference from 42 to 58 cm
Maximum usage time:	15 minutes
Hydrogel composition:	20% water + 30% polyacrylamide + 50% glycerine
Dimensions approx.:	22 × 15 cm
Weight approx.:	520 g
Storage conditions:	0°C – 30°C (32°F – 86°F), protected from light and moisture
Item / EAN numbers:	88110 / 40 15588 88110 7

Disposal

Discard if there are signs of damage, cracks or holes. This product and its packaging must be disposed of in accordance with local regulations.

Warranty and repair terms

Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

- The warranty period for **medisana** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
- Defects in material or workmanship will be removed free of charge within the warranty period.
- Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
- The following is excluded under the warranty:
 - All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - Accessories which are subject to normal wear and tear.
- Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

The current version of this instruction manual can be found under

<https://docs.medisana.com/88110>.

Information about service can be found here:

<https://www.medisana.com/servicepartners>


medisana GmbH,
Carl-Schurz-Str. 2, 41460 NEUSS, GERMANY

Wuhan Pakcare Bio-Technology Co., Ltd.

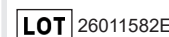Address: No. 28-1, Linkong North Road,
Hengdian Street, Huangpi district, Wuhan city,
430301 Hubei, P.R. China
E-mail: info@pakcare.com

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Share Info GmbH